

ReSound Smart Fit Quick Reference Guide

This guide will help you locate important features as you begin to fit patients. For more detailed instructions, please refer to the ReSound LiNX 3D Fitting Guide

EXPANDED SIDEBAR

Labels on the left side of the image:

- Data logging access
- Connect hearing instruments
- Shows view of hearing instruments associated with this patient
- Mute hearing instruments
- Media Player
- Collapse sidebar

The screenshot shows the 'EXPANDED SIDEBAR' on the left side of the software interface. The sidebar contains the following elements:

- Data logging access:** A button labeled 'Data Logging' at the top of the sidebar.
- Connect hearing instruments:** Two buttons labeled 'Right' and 'Left' with 'Connect' text below them.
- Shows view of hearing instruments associated with this patient:** A visual representation of two hearing aids, one for the right ear and one for the left ear.
- Mute hearing instruments:** A button with a speaker icon and a slash through it, used to mute the hearing instruments.
- Media Player:** A section labeled 'Media Player' with a play button and a volume control slider.
- Collapse sidebar:** A button at the bottom of the sidebar labeled 'Collapse Sidebar'.

The main interface shows the 'Patient' screen with tabs for 'Patient', 'Instruments', 'Fitting', and 'Summary'. The 'Patient' tab is active, displaying 'Patient Experience Level', 'Initial Programs', and 'Right Instrument' and 'Left Instrument' details. An audiogram graph is visible in the top right corner.

PATIENT SCREEN: Profile

Labels on the left side of the image:

- Patient Experience Level
- Initial Programs
- Fit Patient Button

The screenshot shows the 'PATIENT SCREEN: Profile' view. The sidebar is collapsed. The main interface displays the following elements:

- Patient Experience Level:** A dropdown menu showing 'First time user'.
- Initial Programs:** A list of programs including 'All Around', 'Restaurant', 'None', and 'None'.
- Fit Patient Button:** A red button labeled 'Fit Patient' located at the bottom right of the interface.

The main interface also shows 'Right Instrument' and 'Left Instrument' details, including 'National LiNX 3D R' and 'National LiNX 3D L' models, and an audiogram graph.

PATIENT SCREEN: ReSound Assist

Online Services Login

Turn on Remote Fine-tuning or Remote Hearing Aid Updates

FITTING SCREEN: Gain Adjustments

Mute

Display Options

Link/Unlink Hearing Instruments

Gain Level Adjustment

Program Tabs

Gain Adjustment Area

Save

FITTING SCREEN: Advanced Features

Tools Include:

- Reorder Programs
- Recalculate
- AutoRelate
- Manage Program Names
- Copy Program
- Paste Program

Slider to adjust Advanced Features

Advanced Features

Additional Advanced Features


Save

FITTING SCREEN: Environmental Optimizer




Slider to decrease or increase gain for specific environments

FITTING SCREEN: Device Controls

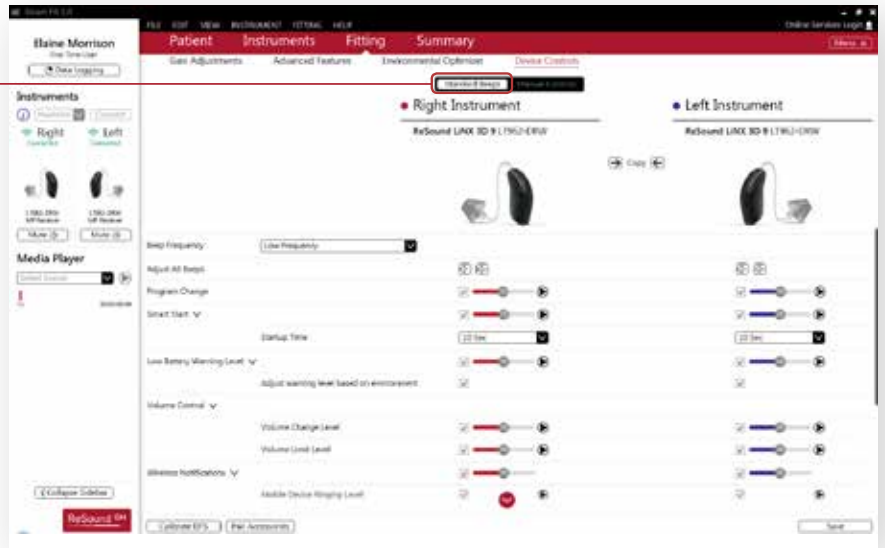

Standard Beeps

FITTING SCREEN: Device Controls

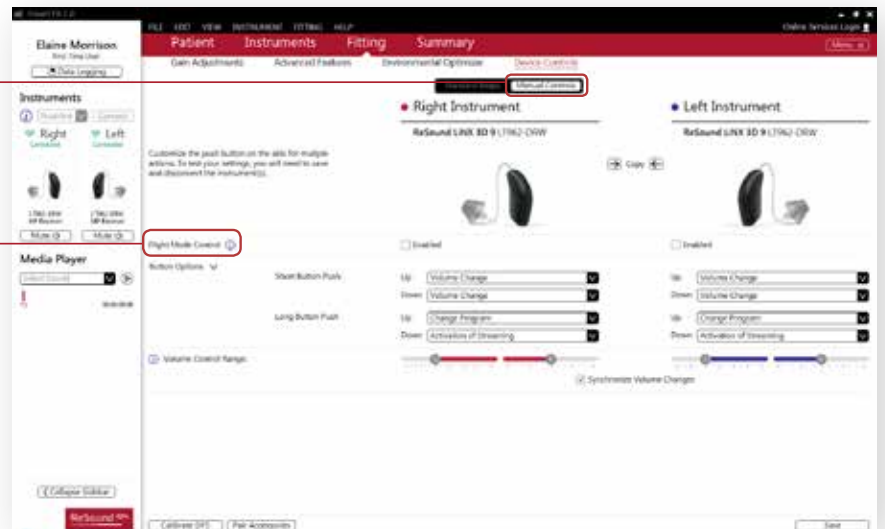


Manual Controls

Enable/Disable Flight Mode


HAMBURGER MENU

Menu Hamburger menu
Houses functions that may not be used during every session

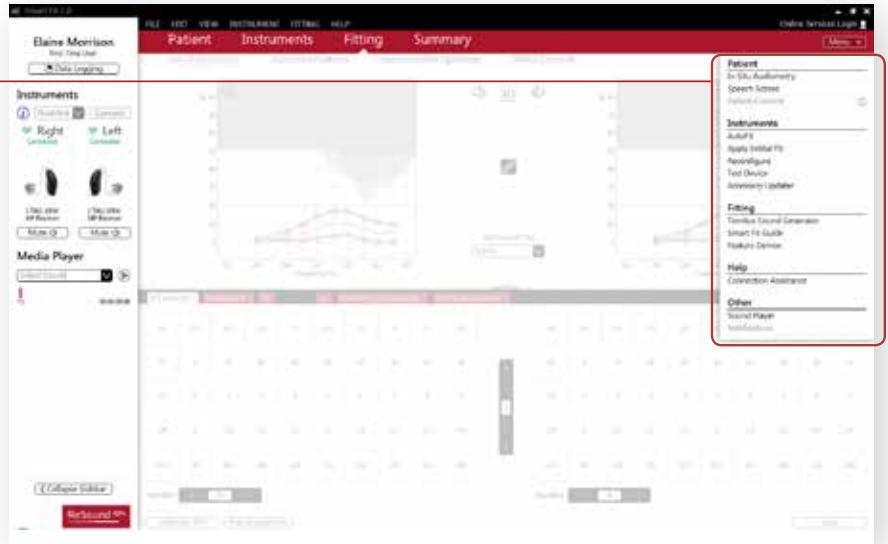

COLLAPSED SIDEBAR

Head icon – Hover to see patient name and patient experience level
Pie chart – Click to open data logging
Red dot by hearing instruments – Not connected
Green dot by hearing instruments – Connected
Play icon – Access media player
Chevron – Expand sidebar

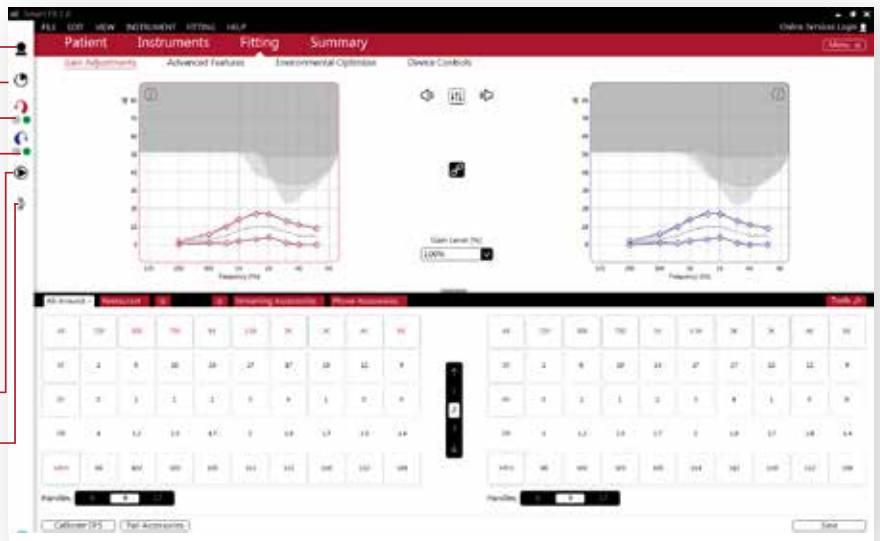

Manufacturer according to FDA:

RESOUND NORTH AMERICA

8001 E. Bloomington Freeway
Bloomington, MN 55420
1.800.248.4327
resoundpro.com

RESOUND GOVERNMENT SERVICES

8001 E. Bloomington Freeway
Bloomington, MN 55420
1.800.392.9932
resound.com/governmentservices

Manufacturer according to Health Canada:

RESOUND CANADA

303 Supertest Road
Toronto, Ontario M3J 2M4
1.888.737.6863
resoundpro.com

