



ReSound GN

Quick guide

ReSound LiNX Quattro™ 67

GN Making Life Sound Better

My hearing aid information

Product: ReSound LiNX Quattro

Model: 67

Serial Number (R):

Serial Number (L):

Battery Size: 312

Push Button:

(if set up by your hearing care professional)

ACTION	RESULT
Short Push	
Long Push	


Programs:

(if set up by your hearing care professional)

PROGRAM	WHEN TO USE	BEEPS
1		1
2		2
3		3
4		4

Get to know your hearing aid

- A. Program Button
- B. Microphone openings
- C. Programming port for cabled fitting
- D. Battery door
- E. Right/Left indicator
- F. Thin tube
- G. Dome

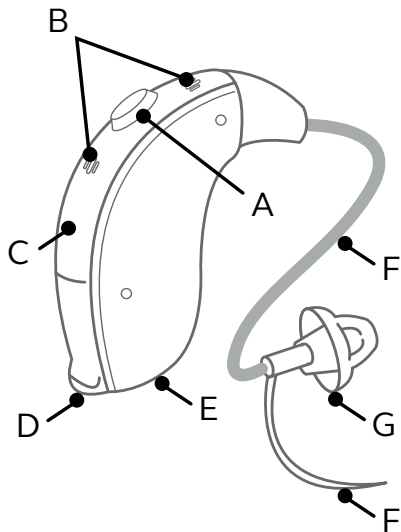


View ReSound hearing aid tutorial videos at youtube.com/resoundus


What you need to know

Getting Started: On/Off function

ON: Insert battery & completely close battery door (D). Program 1 will be active.

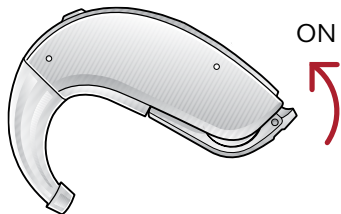

OFF: Open the battery door (D). Use fingertip to pull it all the way.

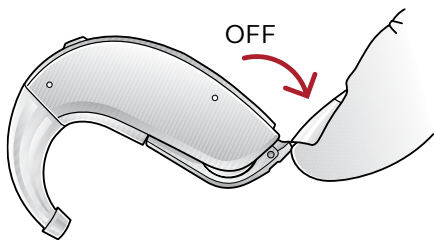

Inserting/Replacing the battery:

1. Open the battery door (D) completely by using your fingertip.



2. Remove protective sticker from battery and wait two minutes before inserting the battery into the hearing aid to allow full activation of the battery.




3. Insert the new battery with the positive side (flat side) in the correct position, as shown above.

4. Gently close the battery door. If it doesn't close easily, check to see if the battery is upside down.

Regular care & maintenance

Daily: Clean your hearing aid.

Keep your hearing aids dry. Clean the hearing aids using a soft cloth or tissue.

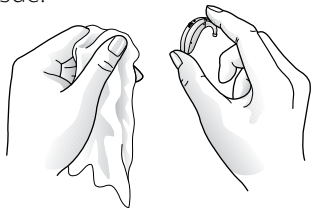


Weekly: Clean the thin tubes and domes/earmold.

1. Remove thin tubes (F) from hearing aids before cleaning by unscrewing them counter clockwise.
2. Wipe down thin tubes and domes (G) with a soft cloth.
3. In order to clear the thin tube of moisture and debris, push the black cleaning rod through the thin tube, beginning at the end opposite the dome.

As needed: Change dome (G) (if applicable)

ReSound standard domes


1. Push the new dome on to the thin tube.
2. Make sure that the new dome is properly and securely maintained.

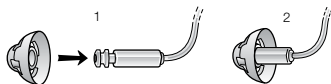

ReSound tulip domes

Tulip domes consist of two “petals”. It is important that the largest petal is the outermost petal. To ensure this:

1. Push the largest petal away from the thin tube using a finger. This bends the petal forward.
2. Then push the largest “petal” back, and it will be placed on top of the smaller petal.

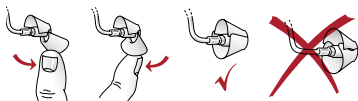

Beep/tone indicators

BEEP/TONE	INDICATION	ACTION REQUIRED
Ding, ding, ding, dong	Low battery	Change battery
Ding (5 or 10)	Hearing aid is starting up	After dings cease, hearing aid is turned on
Ding-ding, ding-ding, etc.	Hearing aid is in flight mode	Use in flight mode if desired, or open-close battery door to exit flight mode
Ripple Tone	Hearing aid is in streaming mode	Use with wireless accessory or press multi-function button to exit streaming mode

View ReSound hearing aid tutorial videos at youtube.com/resoundus

Troubleshooting

SYMPTOM	CAUSE	POSSIBLE REMEDY
No Sound	Battery door not closed	Turn on by closing the battery door (D)
	Dead battery	Replace battery
	Battery door will not close	Insert battery properly
	Blocked earmold or tube	Clean earmold or tube (F)
	Blocked wax filter	Replace wax filter or consult your hearing care professional
Not loud enough	Blocked earmold or dome	Clean earmold, replace dome (G)
	Change in hearing, volume too low	Consult your hearing care professional
	Excessive ear wax	Consult your hearing care professional
Excessive whistling/feedback	Incorrect earmold placement	Re-insert earmold carefully
	Incorrect dome placement	Re-insert dome (G)
	Excessive ear wax	Clean guard/consult your hearing care professional
	If none of the above	Consult your hearing care professional
Wireless does not work	Device may be in flight mode	Open and close the battery door (D) once. Wireless operation will reactivate ten seconds later.

If there are any other problems not mentioned in this guide, please contact ReSound or your hearing care professional.

Flight Mode:

Some airlines prohibit the use of wireless hearing aids (RF transmitters) onboard. If required by your airline, enable flight mode:

Disable wireless operation by opening and closing the battery door (D) three times within ten seconds (open-close, open-close, open-close). Your instruments will now be in flight mode.

The hearing aids must function in flight mode for ten seconds before wireless operation can be re-enabled.

To re-enable wireless operation, open and close the battery door one time and wait thirty seconds before opening the battery door again to avoid accidentally entering flight mode again.

For more information

Refer to the user guide included in your hearing aid box.

Your hearing care professional can be reached at:

View ReSound hearing aid tutorial videos at youtube.com/resoundus.

Manufacturer according to FDA:

ReSound North America

8001 E. Bloomington Freeway
Bloomington, MN 55420
USA
1-888-735-4327
resound.com

ReSound Government Services

8001 E. Bloomington Freeway
Bloomington, MN 55420
USA
1-800-392-9932
www.resound.com/veterans



401298011

401298011 Rev A 2020.01