

## GETTING STARTED *(Quick guide)*

### Phone Clip+

Before you can use your Phone Clip+ it is necessary to:

- A. Charge your Phone Clip+
- B. Pair it with your hearing instrument
- C. Pair it with your mobile phone
- D. Connect it to your mobile phone

Consult your hearing healthcare professional to determine which, if any, of these procedures may have already been completed for you.

## GETTING STARTED *(Quick Guide)*

### Phone Clip+



### **A** CHARGE YOUR PHONE CLIP+

1. Plug one end of the charging cable into the Phone Clip+
2. Plug the other end of the charging cable into an electric wall outlet.



3. Charging is complete when the red LED indicator at the top of the Phone Clip+ turns green with a red blink every 2 seconds. When charging for the first time, the Phone Clip+ should be charged for at least 3 hours, even if the LED turns green before 3 hours have elapsed.



### **B** PAIR WITH YOUR HEARING INSTRUMENT

1. Turn ON your Phone Clip+ by sliding the ON/OFF switch down to the ON position. A white "I" icon will be visible on the base of the switch and the LED indicator will now blink green every 2.5 seconds.
2. Open the battery door on your hearing instruments




### **B** *Pair with your hearing instrument... continued*

3. Remove the silver cap from the Phone Clip+.
4. Press the small white pairing button to put the Phone Clip+ into pairing mode for 20 seconds. While in pairing mode the LED will blink yellow every 2 seconds.



5. While in pairing mode, turn on the hearing instruments by closing the battery doors on both hearing instruments.
6. A melody of beeps will be heard in both hearing instruments once they successfully pair with the Phone Clip+. Replace the cap on the Phone Clip+ when pairing is complete.

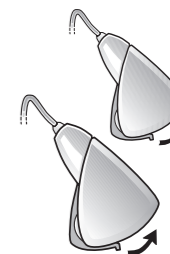

400231244US-17.08-Rev.B

Manufacturer according to FDA:

ReSound North America  
8001 Bloomington Freeway  
Bloomington, MN 55420  
1-888-735-4327

Manufacturer according to Health Canada:

GN Hearing A/S  
Lautrupbjerg 7  
DK-2750 Ballerup

Distributor in Canada:

ReSound Canada  
303 Supertest Road  
Toronto, Ontario  
M3J 2M4  
1-888-737-6863

## C PAIR WITH YOUR MOBILE PHONE

1. Turn your Phone Clip+ ON and remove the silver cap.
2. Find the connectivity settings on your phone. Look for the “Bluetooth” feature in the main menu or one of the submenus, such as “Connectivity”. (Refer to your mobile phone user guide or [www.resound.com/PairPhone](http://www.resound.com/PairPhone) for more detailed instructions.)
3. On your mobile phone, ensure that Bluetooth is turned ON.
4. Press the small blue Bluetooth pairing button one time to activate Bluetooth pairing mode for two minutes. The LED indicator will turn blue during Bluetooth pairing.




5. On the mobile phone search for Bluetooth devices.
6. Your mobile phone should present a list of Bluetooth devices discovered. Select “Hearing Aid Phone” from this list.
7. Your mobile phone may then ask you to enter a passcode. If so, enter “0000” (four zeros).

Your Phone Clip+ should now be paired to your mobile phone. Some mobile phones may ask which Bluetooth service you wish to enable. Select “Headset”.

## D CONNECT WITH YOUR MOBILE PHONE








1. In most instances, the phone will ask you to confirm your intention to connect with the Phone Clip+ after successfully pairing.
2. If you have to manually connect: Locate a list of “Paired Devices” in your phone’s Bluetooth” menu. From that list, highlight “Hearing Aid Phone” and select “Connect.”
3. You are now ready to use your Phone Clip+.
4. On some phones you can configure “Automatic Reconnection.” Refer to the user guide of your mobile phone for instructions on how to do this.

**Tip:** In most instances, a headset symbol  on the main screen of your mobile phone will confirm that you have successfully connected.

**Note:** This Quick Guide is only valid together with the full user guide delivered with your Phone Clip+.

## E LIGHT INDICATORS (LED)

The Light indicator at the top of your Phone Clip+ serves as a multipurpose user interface, providing information on the status of your device.

Blinks	Meaning
	Normal operation
	Battery low
	Charging
	Fully charged -still connected to charger
	Ready for pairing with hearing instrument(s)
	Ready for pairing with Bluetooth device
	Bluetooth call active

## F PHONE OPERATION




In the below table you will find an overview of how you operate the Phone Clip+ when using it in connection with phone calls.

<b>Answering and ending a call</b>	Push the call pick-up/hang-up button once
<b>Rejecting a call</b>	Double click the call pick-up/hang-up button
<b>Transferring a call initiated on the mobile phone to Phone Clip+</b>	Push the call pick-up/hang-up button once
<b>Put an active call on hold</b>	Press and hold the call pick-up/hang-up button for 2 seconds
<b>Last number redial</b>	Double click the call pick-up/hang-up button while no call is active nor pending
<b>Voice dialing</b>	Press and hold the call pick-up/hang-up button for 2 seconds

**Note:** Except for the answer and end call operations, one or more of the below described operations may be phone dependent, i.e. they may either differ slightly or be absent altogether.

## G REMOTE CONTROL OPERATION

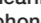When not on the phone, the Phone Clip+ can be used as a simple remote control for your hearing instruments\*.

<b>Turn hearing instrument(s) volume up</b>	Push + on the volume button
<b>Turn hearing instrument(s) volume down</b>	Push - on the volume button
<b>Mute hearing instrument microphones</b>	Push  once
<b>Unmute hearing instrument microphones</b>	Push  once
<b>Toggle through hearing instrument programs</b>	Push  one or more times

\*Unless upgraded to the latest software, some wireless hearing instruments must be initially paired with a compatible Remote Control to enable the Phone Clip’s remote control feature. This may be done with either your or your hearing healthcare professional’s Remote Control. Mute and stereo streaming is unavailable on certain hearing instrument models. Consult your hearing healthcare professional for more information.