

# Charging options for Beltone hearing aids

Beltone offers rechargeable solutions with inductive charging and long battery life to keep your clients on the go.



Premium charger







Standard charger



Desktop charger



Custom charger

Beltone Achieve™ miniRIE		06	06	06	-
Beltone Achieve™ BTE Beltone Imagine™ BTE		-	-	70	-
Beltone Achieve™ PBTE Beltone Imagine™ PBTE		-	-	80	-
Beltone Achieve™ ITC/ITE Beltone Imagine™ ITC/ITE		-	-	-	C

## How to match Beltone hearing aids with the correct charger

Each charger type is offered with different inserts – known as the charging form factor (CFF) – which relate to a specific Beltone hearing aid collection or style.

The CFF is indicated with a number and is shown inside the charger and on the charger packaging.



**CFF numbers are not shown on all charger packaging yet, please refer the number inside the charger:**

- Beltone Imagine RIE & Beltone Achieve RIE have CFF 1 shown inside the charger
- Beltone Amaze RIE & Beltone Rely RIE have CFF 4 shown inside the charger



# Charger features

	Beltone Achieve			Beltone Achieve Beltone Imagine						Beltone Amaze Beltone Rely	
	Premium charger	Standard charger	Desktop charger	Premium charger	Standard charger	Desktop charger	Desktop charger	Desktop charger	Custom charger	Premium charger	Desktop charger
	miniRIE	miniRIE	miniRIE	RIE	RIE	RIE	BTE	Power BTE	Customs	RIE	RIE
Inductive charging	•	•	•	•	•	•	•	•	•	•	•
Portable with onboard power bank	•			•						•	
Portable without power bank		•			•				•		
USB-C connector and cable	•	•	•	•	•	•	•	•	•	•	•
LED on the back	Shows the power bank battery status	Lights up when the charger is connected to the main power supply	Shows when the charger is connected to a power supply	Shows the power bank battery status	Lights up when the charger is connected to the main power supply	Shows when the charger is connected to a power supply	Shows when the charger is connected to a power supply	Shows when the charger is connected to a power supply	Shows when the charger is connected to a power supply	Shows the power bank battery status	Shows when the charger is connected to a power supply
LED on the front	Shows the hearing aids charging			Shows the hearing aids charging					Shows the hearing aids charging	Shows the hearing aids charging	

## How to tell if the hearing aids are charging

Place the hearing aids in the charger and let the lights indicator on the hearing aids guide you. It is the same process whichever charger is used.

You know that charging has started when the hearing aid lights blink slowly.

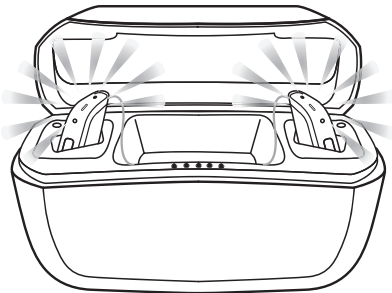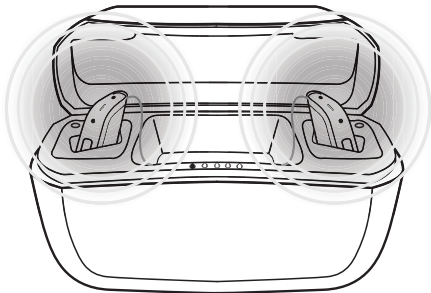

When fully charged, the hearing aid lights remain constant, not blinking.

The lights turn off when the hearing aids are removed from the charger.

Premium charger illustrated.  
For more information, please see the individual user guide for each charger.

Worldwide Headquarters  
Beltone A/S  
Lautrupbjerg 7  
DK-2750 Ballerup, Denmark  
Tel.: +45 45 75 11 11  
beltone.com

CVR no. 55082715

Australia  
GN Hearing Australia Pty Ltd  
Gate C, 19-25 Khartoum Road  
Macquarie Technology Park  
Macquarie Park NSW 2113  
Tel.: (free) 1800 658 955  
beltone.com/en-au

New Zealand  
GN Hearing New Zealand Limited Ground  
Ground Floor, North Entrance  
4 Fred Thomas Drive  
Takapuna, Auckland, 0622  
Tel.: (free) 0800 900 126  
beltone.com/en-au

