

# Charging options for ReSound hearing aids

ReSound offers rechargeable solutions with inductive charging and all-day battery life to keep your clients on the go.








Premium charger

Standard charger

Desktop charger

Custom charger

ReSound Nexia™ microRIE 	8	8	8	-
ReSound OMNIA™ miniRIE 	06	06	06	-
ReSound OMNIA™ BTE 	-	-	70	-
ReSound OMNIA™ PBTE 	-	-	80	-
ReSound OMNIA™ ITE/ITC 	-	-	-	C

## How to match ReSound hearing aids with the correct charger

Each charger type is offered with different inserts – known as the charging form factor (CFF) – which relate to a specific ReSound hearing aid collection or style.

The CFF is indicated with a number and is shown inside the charger and on the charger packaging.



**CFF numbers are not shown on all charger packaging yet, please refer the number inside the charger:**

- ReSound OMNIA RIE have CFF 1 shown inside the charger
- ReSound Quattro RIE & ReSound Key RIE have CFF 3 shown inside the charger



# Charger features

	ReSound Nexia			ReSound OMNIA				ReSound LINX Quattro ReSound Key	
	Premium charger	Standard charger	Desktop charger	Premium charger	Standard charger	Desktop charger	Custom charger	Premium charger	Desktop charger
	microRIE	microRIE	microRIE	miniRIE/RIE	miniRIE/RIE	miniRIE/RIE BTE/Power BTE	Customs	RIE	RIE
Inductive charging	•	•	•	•	•	•	•	•	•
Portable with onboard power bank	•			•				•	
Portable without power bank		•			•		•		
USB-C connector and cable	•	•	•	•	•	•	•	•	•
LED on the back	Shows the power bank battery status	Lights up when the charger is connected to the main power supply	Shows when the charger is connected to a power supply	Shows the power bank battery status	Lights up when the charger is connected to the main power supply	Shows when the charger is connected to a power supply	Shows when the charger is connected to a power supply	Shows the power bank battery status	Shows when the charger is connected to a power supply
LED on the front	Shows the hearing aids charging			Shows the hearing aids charging			Shows the hearing aids charging	Shows the hearing aids charging	

## How to tell if the hearing aids are charging

Place the hearing aids in the charger and let the lights indicator on the hearing aids guide you. It is the same process whichever charger is used.

You know that charging has started when the hearing aid lights blink slowly.

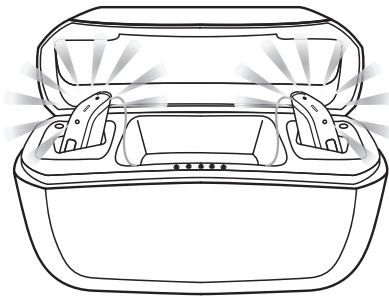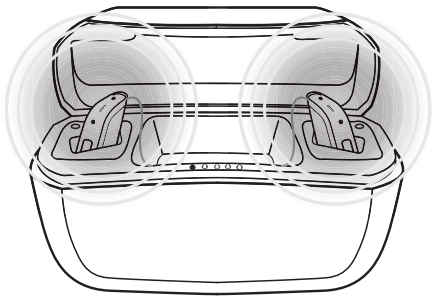

When fully charged, the hearing aid lights remain constant, not blinking.

The lights turn off when the hearing aids are removed from the charger.

Premium charger illustrated.  
For more information, please see the individual user guide for each charger.

**Worldwide Headquarters**  
GN ReSound A/S  
Lautrupbjerg 7  
DK-2750 Ballerup  
Denmark  
Tel.: +45 4575 1111  
resound.com  
CVR no. 55082715

**United Kingdom**  
GN Hearing UK Ltd.  
Unit 13  
Talisman Business Centre  
Bicester  
Oxon OX26 6HR  
Tel: +44 1869 352800  
resound.com

**Australia**  
GN Hearing Australia Pty Ltd  
Gate C, 19-25 Khartoum Road  
Macquarie Technology Park  
Macquarie Park NSW 2113  
Tel.: (free) 1800 658 955  
resound.com

**New Zealand**  
GN Hearing New Zealand Limited  
Ground Floor, North Entrance  
4 Fred Thomas Drive  
Takapuna, Auckland, 0622  
Tel.: (free) 0800 900 126  
resound.com